

Thermo Scientific ClipTip Pipette Tips

Thermo Scientific™ ClipTip™ pipette tips provide a unique locking tip interface that allows stable, straight mounting on Thermo Scientific F1-ClipTip manual pipettes or E1-ClipTip electronic pipettes. ClipTip pipette tips feature “clips” that lock tips in place with a light touch for a secure seal on every channel. ClipTip pipette tips will only release when ejected using minimal application pressure.

Absolute tip attachment.

The ClipTip pipette tips lock into place to create a complete seal until released. Users have full tip control throughout pipetting, as pipette tips will only release when ejected, regardless of application pressure.

‘Snap and Seal’ tip attachment technology.

The 96-format ClipTip 12.5 and all ClipTip 384 tips feature a ‘Snap and Seal’ tip attachment mechanism. The tip has incorporated small protrusions that snap on the tip fitting flange, ensuring light-force and secure tip attachment.

Superior ergonomics.

You only need a light touch to attach ClipTip tips on the F1-ClipTip or E1-ClipTip pipettes to lock them in place.

Immediate traceability.

Each tip insert has ink-jetted lot and volume information. The convenient color key makes it easy to match the color of the F1-ClipTip pipette plunger cap or buttons on the E1-ClipTip pipette with the corresponding ClipTip pipette tip rack insert.



Superior Contamination Control.

The interlocking tip design prevents sample contamination caused from loose, leaking tips that can fall off. Filter tips prevent carry over contamination protecting both the pipette and the sample.

Low Retention.

Innovative technology creates a hydrophobic inner pipette tip surface to minimize liquid retention, and enable maximum sample recovery. Available in ClipTip 20-1250 96-format tips.

Proven Quality.


Certified free of RNase, DNase, DNA, ATP and Endotoxin.

Environmental responsibility.

The ClipTip pipette tip reload system incorporates re-usable hinged racks with plastic reload sleeves that reduces waste by an average of 35%.

Versatile Tip Range.

- Choose from seven sizes in 96-format: 12.5, 20, 50, 200, 300, 1000 and 1250 μ L
- Choose from three sizes in 384-format: 12.5, 30 and 125 μ L

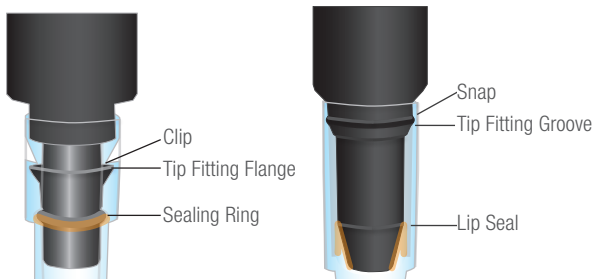

ClipTip interlocking technology that locks tips into place for a complete seal every time.

The 96-format ClipTip 12.5 and all ClipTip 384 tips feature a ‘Snap and Lip’ seal mechanism for tip attachment.

Thermo Scientific ClipTip Pipette Tips


Item No.	Description	Color Key	Volume Range	Packaging
96-format ClipTip Tips				
94410040	ClipTip 12.5*	●	0.5-12.5 µL	10 x 96/rack
94410043	ClipTip 12.5, Sterile*	●	0.5-12.5 µL	10 x 96/rack
94420043	ClipTip 12.5, Filter, Sterile*	●	0.5-12.5 µL	10 x 96/rack
94410210	ClipTip 20	●	1-20 µL	10 x 96/rack
94410213	ClipTip 20, Sterile	●	1-20 µL	10 x 96/rack
94410217	ClipTip 20, Reload Tower	●	1-20 µL	10 x 96/insert
94410218	ClipTip 20, Reload Tower, Sterile	●	1-20 µ	10 x 96/insert
94420213	ClipTip 20, Filter, Sterile	●	1-20 µL	10 x 96/rack
94420218	ClipTip 20, Filter Reload Insert, Sterile	●	1-20 µL	10 x 96/insert
94410250	ClipTip 50	●	5-50 µL	10 x 96/rack
94410253	ClipTip 50, Sterile	●	5-50 µL	10 x 96/rack
94410257	ClipTip 50, Reload Tower	●	5-50 µL	10 x 96/insert
94410258	ClipTip 50, Reload Tower, Sterile	●	5-50 µL	10 x 96/insert
94420253	ClipTip 50, Filter, Sterile	●	5-50 µL	10 x 96/rack
94420258	ClipTip 50, Filter Reload Insert, Sterile	●	5-50 µL	10 x 96/insert
94410310	ClipTip 200	●	2-200 µL	10 x 96/rack
94410313	ClipTip 200, Sterile	●	2-200 µL	10 x 96/rack
94410317	ClipTip 200, Reload Tower	●	2-200 µL	10 x 96/insert
94410318	ClipTip 200, Reload Tower, Sterile	●	2-200 µL	10 x 96/insert
94420313	ClipTip 200, Filter, Sterile	●	2-200 µL	10 x 96/rack
94420318	ClipTip 200, Filter Reload Insert, Sterile	●	2-200 µL	10 x 96/insert
94410510	ClipTip 300	●	0-300 µL	10 x 96/rack
94410513	ClipTip 300, Sterile	●	10-300 µL	10 x 96/rack
94410517	ClipTip 300, Reload Tower	●	10-300 µL	10 x 96/insert
94410518	ClipTip 300, Reload Tower, Sterile	●	10-300 µL	10 x 96/insert
94420513	ClipTip 300, Filter, Sterile	●	10-300 µL	10 x 96/rack
94420518	ClipTip 300, Filter Reload Insert, Sterile	●	10-300 µL	10 x 96/insert
94410710	ClipTip 1000	●	30-1000 µL	8 x 96/rack
94410713	ClipTip 1000, Sterile	●	30-1000 µL	8 x 96/rack
94410717	ClipTip 1000, Reload Tower	●	30-1000 µL	8 x 96/insert
94410718	ClipTip 1000, Reload Tower, Sterile	●	30-1000 µL	8 x 96/insert
94420713	ClipTip 1000, Filter, Sterile	●	30-1000 µL	8 x 96/rack
94420718	ClipTip 1000, Filter Reload Insert, Sterile	●	30-1000 µL	8 x 96/insert
94410810	ClipTip 1250	●	15-1250 µL	8 x 96/rack
94410813	ClipTip 1250, Sterile	●	15-1250 µL	8 x 96/rack
94410817	ClipTip 1250, Reload Insert	●	15-1250 µL	8 x 96/insert
94410818	ClipTip 1250, Reload Insert, Sterile	●	15-1250 µL	8 x 96/insert
94420813	ClipTip 1250, Filter, Sterile	●	15-1250 µL	8 x 96/rack
94420818	ClipTip 1250, Filter Reload Insert, Sterile	●	15-1250 µL	8 x 96/insert
384 ClipTip Tips*				
94410050	ClipTip 384 12.5	●	0.5-12.5 µL	10x384/rack
94410053	ClipTip 384 12.5, Sterile	●	0.5-12.5 µL	10x384/rack
94420053	ClipTip 384 12.5, Filter, Sterile	●	0.5-12.5 µL	10x384/rack
94410100	ClipTip 384 30	●	1-30 µL	10x384/rack
94410103	ClipTip 384 30, Sterile	●	1-30 µL	10x384/rack
94420103	ClipTip 384 30, Filter, Sterile	●	1-30 µL	10x384/rack
94410150	ClipTip 384 125	●	2-125 µL	10x384/rack
94410153	ClipTip 384 125, Sterile	●	2-125 µL	10x384/rack
94420153	ClipTip 384 125, Filter, Sterile	●	2-125 µL	10x384/rack

* ClipTip 12.5 and ClipTip 384 tips do not have interlocking tip interface or low retention technology

Thermo Scientific ClipTip Pipette Tip Dimensions

Description	Overall Tip Length	Tip Shoulder (Proximal)	Inner Orifice Diameter (Distal)
96 ClipTip Tips			
ClipTip 12.5*	30 mm	9 mm	0.4 mm
ClipTip 20	46 mm	17 mm	0.4 mm
ClipTip 50	51 mm	17 mm	0.4 mm
ClipTip 200	56 mm	22 mm	0.5 mm
ClipTip 300	63 mm	22 mm	0.6 mm
ClipTip 1000	95 mm	32 mm	1.1 mm
ClipTip 1250	103 mm	32 mm	0.9 mm
384 ClipTip Tips*			
ClipTip 384 12.5	30 mm	9 mm	0.4 mm
ClipTip 384 30	44 mm	9 mm	0.4 mm
ClipTip 384 125	64 mm	9 mm	0.5 mm

* All ClipTip pipette tips are low retention with the exception of 384 and 12.5 models. ClipTip 384 and ClipTip 12.5 tips and compatible E1-ClipTip models do not have interlock technology.

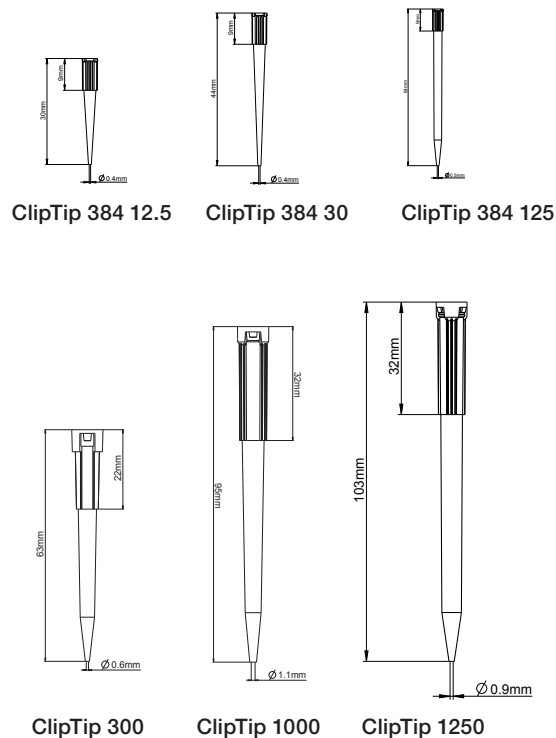

ClipTip Compatibility

Order No.	Description						
		ClipTip 20	ClipTip 50	ClipTip 200	ClipTip 300	ClipTip 1000	
F1-ClipTip Manual Pipettes							
Single channels							
4641170N	F1-ClipTip 1-10 µL	•					
4641180N	F1-ClipTip 2-20 µL	•					
4641190N	F1-ClipTip 5-50 µL		•				
4641200N	F1-ClipTip 10-100 µL			•			
4641210N	F1-ClipTip 20-200 µL			•			
4641220N	F1-ClipTip 30-300 µL				•		
4641230N	F1-ClipTip 100-1000 µL					•	
4651190N	F1-ClipTip Fixed 10 µL	•					
4651200N	F1-ClipTip Fixed 20 µL	•					
4651210N	F1-ClipTip Fixed 25 µL		•				
4651220N	F1-ClipTip Fixed 50 µL		•				
4651230N	F1-ClipTip Fixed 100 µL			•			
4651240N	F1-ClipTip Fixed 200 µL			•			
4651250N	F1-ClipTip Fixed 250 µL				•		
4651260N	F1-ClipTip Fixed 500 µL					•	
4651270N	F1-ClipTip Fixed 1000 µL					•	
Multichannels							
4661110N	F1-ClipTip 8-ch 1-10 µL	•					
4661120N	F1-ClipTip 8-ch 5-50 µL		•				
4661130N	F1-ClipTip 8-ch 10-100 µL			•			
4661140N	F1-ClipTip 8-ch 30-300 µL				•		
4661150N	F1-ClipTip 12-ch 1-10 µL	•					
4661160N	F1-ClipTip 12-ch 5-50 µL		•				
4661170N	F1-ClipTip 12-ch 10-100 µL			•			
4661180N	F1-ClipTip 12-ch 30-300 µL				•		

ClipTip Compatibility

Order No.	Description	ClipTip 12.5	ClipTip 20	ClipTip 50	ClipTip 200	ClipTip 300	ClipTip 1000	ClipTip 1250	ClipTip 12.5 384	ClipTip 30 384	ClipTip 125 384
E1-ClipTip Electronic Pipettes											
Single channels											
4670000	E1-ClipTip 0.5-12.5 µL	•									
4670020	E1-ClipTip 2-125 µL				•						
4670030	E1-ClipTip 10-300 µL					•					
4670040	E1-ClipTip 15-1250 µL							•			
Multichannels											
4671000	E1-ClipTip 8-ch 0.5-12.5 µL	•									
4671040	E1-ClipTip 8-ch 2-125 µL				•						
4671070	E1-ClipTip 8-ch 10-300 µL					•					
4671100	E1-ClipTip 8-ch 15-1250 µL							•			
4671010	E1-ClipTip 12-ch 0.5-12.5 µL	•									
4671050	E1-ClipTip 12-ch 2-125 µL				•						
4671080	E1-ClipTip 12-ch 10-300 µL					•					
4671090	E1-ClipTip 12-ch 30-850 µL						•				
4671020	E1-ClipTip 16-ch 0.5-12.5 µL								•		
4671030	E1-ClipTip 16-ch 1-30 µL									•	
4671060	E1-ClipTip 16-ch 2-125 µL										•
Adjustable Tip Spacing Multichannels											
4672090	E1-ClipTip Equalizer 6-ch 15-1250 µL							•			
4672050	E1-ClipTip Equalizer 8-ch 2-125 µL				•						
4672080	E1-ClipTip Equalizer 8-ch 10-300 µL					•					
4672100	E1-ClipTip Equalizer 8-ch 15-1250 µL							•			
4672010	E1-ClipTip Equalizer 384 8-ch 0.5-12.5 µL								•		
4672030	E1-ClipTip Equalizer 384 8-ch 1-30 µL									•	
4672060	E1-ClipTip Equalizer 384 8-ch 2-125 µL										•
4672020	E1-ClipTip Equalizer 384 12-ch 0.5-12.5 µL								•		
4672040	E1-ClipTip Equalizer 384 12-ch 1-30 µL									•	
4672070	E1-ClipTip Equalizer 384 12-ch 2-125 µL										•



Learn more at thermofisher.com/cliptip

© 2016 Thermo Fisher Scientific Inc. All rights reserved.

All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Specifications, terms and pricing are subject to change. Please consult your local sales representative for details.

Australia +61 39757 4300

Austria +43 1 801 40 0

Belgium +32 53 73 42 41

China +86 21 6865 4588 or
+86 10 8419 3588

France +33 2 2803 2180

Germany national Toll-Free 0800 1 536 376

Germany International +49 6184 90 6000

India Toll-Free 1800 22 8374

India +91 22 6716 2200

Italy +39 02 95059 554

Japan +81 45 453 9220

Netherlands +31 76 579 55 55

New Zealand +64 9 980 6700

Nordic/Baltic/GIS Countries

+358 9 329 10200

Russia +7 812 703 42 15

Spain/Portugal +34 93 223 09 18

Switzerland +41 44 454 12 22

UK/Ireland +44 870 609 9203

USA/Canada +1 866 984 3766

Other Asian countries +852 2885 4613

Countries Not Listed +49 6184 90 6000

ThermoFisher
SCIENTIFIC