




Experience seamless PCR sampling, preparation, and storage with this specially curated bundle.

Fisherbrand Mini Centrifuge

This compact centrifuge offers user-friendly features such as tool-free SnapSpin rotors, quick-spin options, a braking function for quick stops, and a weighted rubberized base for stability and noise reduction.

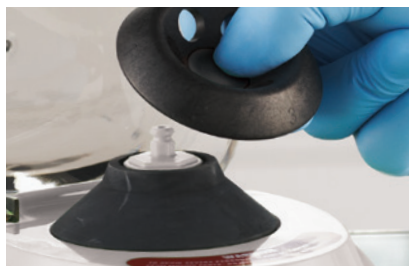

SnapSpin Rotors

change out with the push of a button. No tools are needed!

Fisherbrand Mini Vortexer

Small footprint Vortexer offering maximum vortex shaking for a variety of tube sizes. A simple touch start operation allows for immediate sample mixing with one hand and automatic stop when tube is lifted.

Fisherbrand PCR Cooler

Ice-free cooling system for PCR preparation and storage. Cools industry-standard PCR tubes, strips, and plates with efficiency. Maintains the entire 96-Well plate at a temperature of 0°C for over an hour and keeps samples below 5°C for up to three hours. Vibrant color change signals to user when temperatures exceed 6°C.

Fisherbrand 96 Place Storage Rack Low Temp


These Low Temperature PCR Racks have a rotating sliding lid that doubles as a stand, providing convenience and stability. The translucent lid is specially designed for rapid cooling and quick sample identification. With ergonomic finger grips and a robust ABS base, these racks ensure secure handling and sample safety. Available in multiple colors for lab coding.

Fisherbrand Silicone Lab Mat


This FDA food-grade silicone labmat is suitable for various laboratory settings, including cold rooms, clean rooms, and sterile suites, where high sterility is essential. Effectively safeguards countertops from temperatures up to 200°C (392°F) and demonstrates general chemical resistance to moderate oxidizing chemicals.



Reversible
Blue and Gray

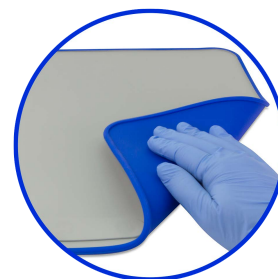

Fisherbrand Mini Centrifuge (16617645)

Colour		L x W x H in	L x W x H cm
Grey		6 x 5 x 4.3	15.3 x 12.8 x 10.9

Specifications			
Speed Range	Approx. 6,000 rpm (fully loaded)		
RCF	2,000 x g		




Package Includes

Mini Centrifuge	6-Place 1.5/2.0 mL Tube Rotor	6 Tube Adapters 0.2 mL
6 Tube Adapters 0.5 mL	4 Lid O-Rings	Tube Rack
Storage Case for Rotors and Adapters		
16-Place 0.2 mL Tube Rotor (singles, or 2 x 8 strip)		
1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs		



Fisherbrand Mini Vortexer (17871672)

Colour		L x W x H in	L x W x H cm
Grey/Blue		5.1 x 4.3 x 3	13 x 11 x 7.5

Specifications	
Speed	4,500 rpm
Orbit	6 mm
Capacity	Up to 50 mL tube



Package Includes


Mini Vortexer	2 Silicone Caps
1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs	



Fisherbrand PCR Cooler (16320412)

Warm Colour	Cold Colour	Wells	Tube Sizes mL	L x W x H in	L x W x H cm
Pink	Purple	96	0.2 Assay Plates	4.3 x 5.8 x 1.9	11 x 14.8 x 4.8
Light Blue	Dark Blue				

Fisherbrand Silicone Lab Mat (15837944)

Colour		L x W x H in	L x W x H cm
Grey/Blue		14 x 23.5 x 0.1	35 x 60 x 0.2

96 Place Storage Rack Low Temp (15758177)

Colour	Tubes mL	Wells	Well Dia. mm	L x W x H in	L x W x H cm
Assorted*	0.2	96	5.2	4.9 x 3.6 x 1.5	12.5 x 9.2 x 3.7

*Assorted colours (Purple, Blue, Green, Orange, Yellow)



Distributed by Fisher Scientific. Contact us today:

Austria: fishersci.at
 Belgium: fishersci.be
 Denmark: fishersci.dk
 Germany: fishersci.de
 Ireland: fishersci.ie
 Italy: fishersci.it
 Finland: fishersci.fi
 France: fishersci.fr
 Netherlands: fishersci.nl
 Norway: fishersci.no
 Portugal: fishersci.pt
 Spain: fishersci.es
 Sweden: fishersci.se
 Switzerland: fishersci.ch
 UK: fishersci.co.uk

© 2023 Thermo Fisher Scientific Inc. All rights reserved.
 Trademarks used are owned as indicated at fishersci.com/trademarks.